



ANALIZA BIAŁEK



100011762

## Laboratory of Analysis

**Client:** HEPiUS SHOP  
**Order number:** 100011762  
**Name:** STROMUSC CARDARINE  
**Number of samples:** 1 container

**Batch number:** -----  
**Receive date:** 2026-01-02  
**Start date:** 2026-01-23  
**End date:** 2026-01-27

### Sample identification

Substance	Method of taking a sample	Form
GW-501516	Sample provided by customer	Blue tablets in a plastic container

### Results

Substance	Sample name	Content	Purity
GW-501516	STROMUSC CARDARINE 10MG/TAB	10,51mg/tab	-----

### Research comment

The test sample was applied in the amount of 3 µl to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H2O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 10,51mg/tab. Height the measurement error is + - 0,05mg. The actual state of the sample sent before the test is presented on the photo below

### Photo of the test sample



### Password

K83NLF9S

**Analiza Białek sp. z o.o.**  
**KRS:** 0000966168  
**NIP:** 8992922970  
**REGON:** 521751875

**Headquarters and Laboratory:**  
ul. Duńska 9, 54-427 Wrocław  
Budynek DELTA

**Lab technician:**  
Analysis and research specialist  
mgr Arkadiusz Zająć




ANALIZA BIAŁEK



100011762

## Chromatogram

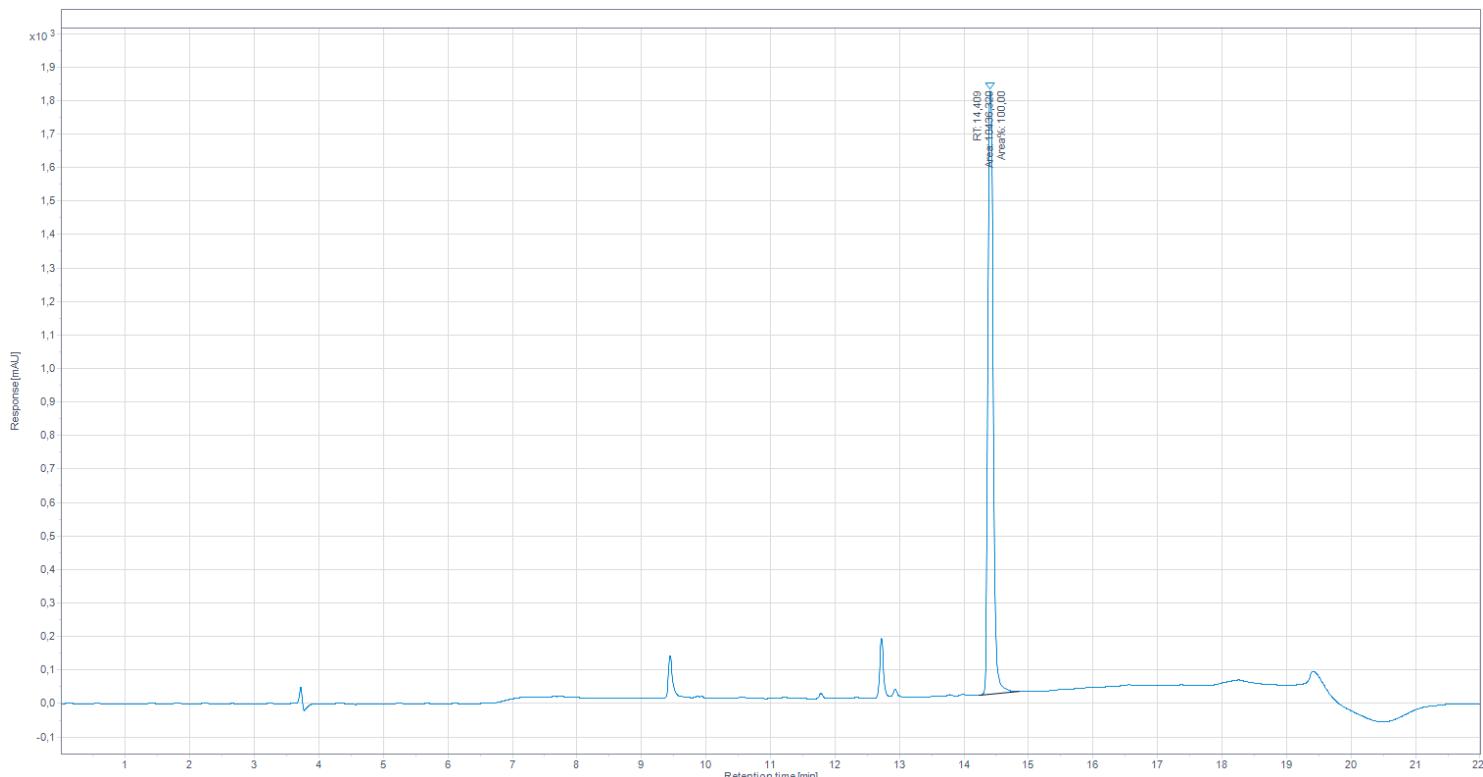

## Method

Time	H <sub>2</sub> O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
12min	5%	95%
15min	5%	95%

## Password

K83NLF9S

**Analiza Białek sp. z o.o.**  
**KRS:** 0000966168  
**NIP:** 8992922970  
**REGON:** 521751875

**Headquarters and Laboratory:**  
ul. Duńska 9, 54-427 Wrocław  
Budynek DELTA

**Lab technician:**  
Analysis and research specialist  
mgr Arkadiusz Zająć