



ANALIZA BIAŁEK



100011768

## Laboratory of Analysis

**Client:** HEPiUS SHOP  
**Order number:** 100011768  
**Name:** STROMUSC STANOZOLOL  
**Number of samples:** 1 container

**Batch number:** -----  
**Recive date:** 2026-01-02  
**Start date:** 2026-01-23  
**End date:** 2026-01-27

### Sample identification

Substance	Method of taking a sample	Form
Stanozolol	Sample provided by customer	White tablets in a plastic container


### Results

Substance	Sample name	Content	Purity
Stanozolol	STROMUSC STANOZOLOL 20MG/TAB	19,57mg/tab	-----

### Research comment

The test sample was applied in the amount of 3 µl to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H<sub>2</sub>O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peak of the sought substance appeared. Its concentration, when compared to the reference chromatogram, is 19,57mg/tab. Height the measurement error is + - 0,05mg. The actual state of the sample sent before the test is presented on the photo below

### Photo of the test sample

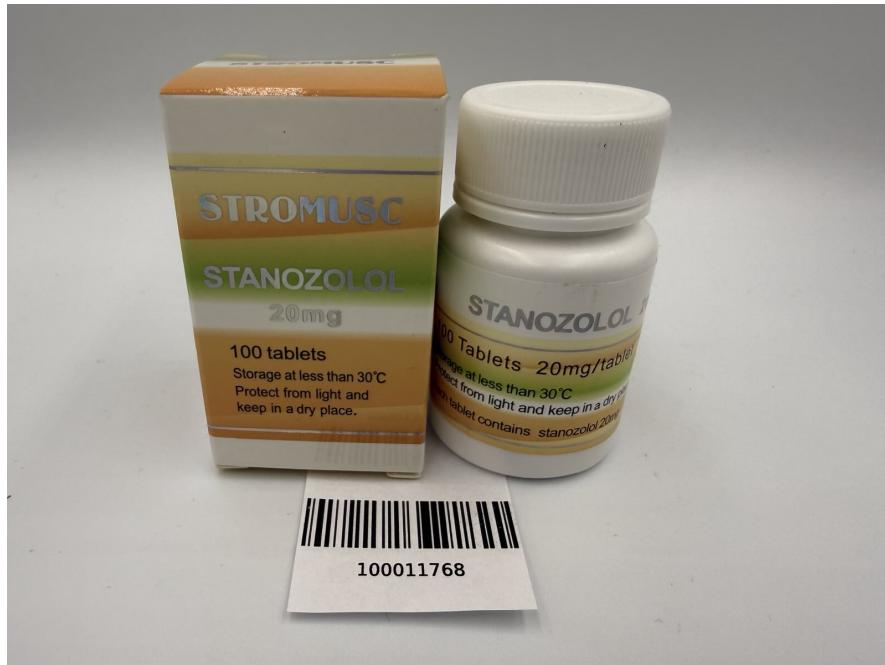

### Password

7RL93NE4

**Analiza Bialek sp. z o.o.**  
**KRS:** 0000966168  
**NIP:** 8992922970  
**REGON:** 521751875

**Headquarters and Laboratory:**  
ul. Duńska 9, 54-427 Wrocław  
Budynek DELTA

**Lab technician:**  
Analysis and research specialist  
mgr Arkadiusz Zająć




ANALIZA BIAŁEK



100011768

## Chromatogram

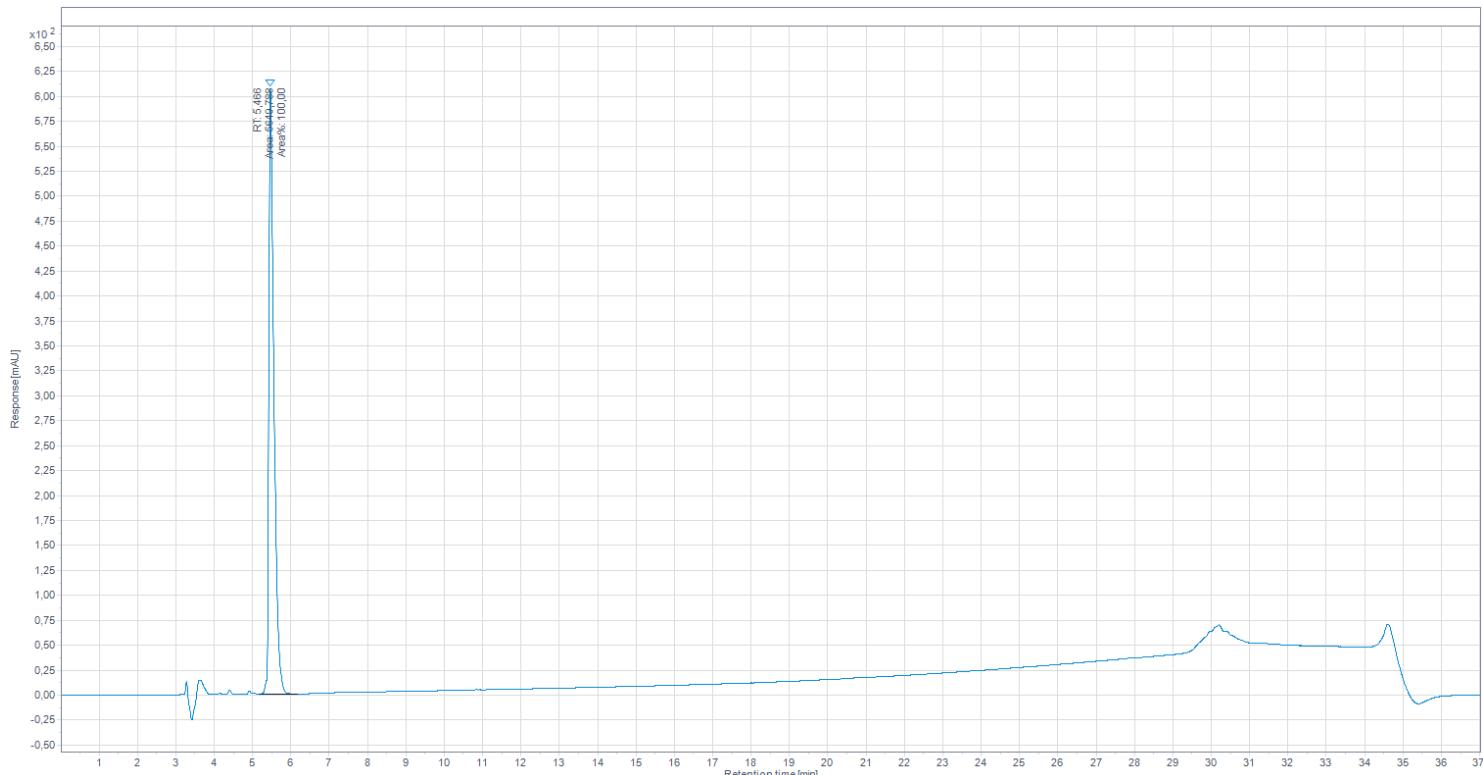

## Method

Time	H <sub>2</sub> O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
12min	5%	95%
15min	5%	95%

## Password

7RL93NE4

**Analiza Białek sp. z o.o.****KRS:** 0000966168**NIP:** 8992922970**REGON:** 521751875**Headquarters and Laboratory:**ul. Duńska 9, 54-427 Wrocław  
Budynek DELTA**Lab technician:**Analysis and research specialist  
mgr Arkadiusz Zająć