

ANALIZA BIAŁEK



100009833

Laboratory of Analysis

Client: HEPIUS SHOP
Order number: 100009833
Name: STROMUSC COCK
Number of samples: 1 container

Batch number: -----
Recive date: 2025-09-17
Start date: 2025-10-09
End date: 2025-10-13

Sample identification

Substance	Method of taking a sample	Form
Tadalafil	Sample provided by customer	White tablets in a plastic container

Results

Substance	Sample name	Content	Purity
Tadalafil	STROMUSC COCK 70MG/TAB	21,60mg/tab	-----
Sildenafil	STROMUSC COCK 70MG/TAB	53,21mg/tab	-----

Research comment

The test sample was applied in the amount of 2 µl to the BioBasic 8 Thermo Scientific column. During the course of the analysis with the method a linear gradient in the solutions of H₂O and Acetonitrile with the addition of 0.1% trifluoroacetic acid, a peaks of the sought substances appeared. Its concentration, when compared to the reference chromatogram shows the table above. Height the measurement error is + - 0,05mcg. The actual state of the sample sent before the test is presented on the photo below

Photo of the test sample




Password

D63NFU74

Analiza Białek sp. z o.o.
KRS: 0000966168
NIP: 8992922970
REGON: 521751875

Headquarters and Laboratory:
ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:
Analysis and research specialist
mgr Arkadiusz Zajac

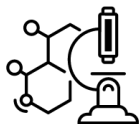


ANALIZA BIAŁEK



100009833

Chromatogram

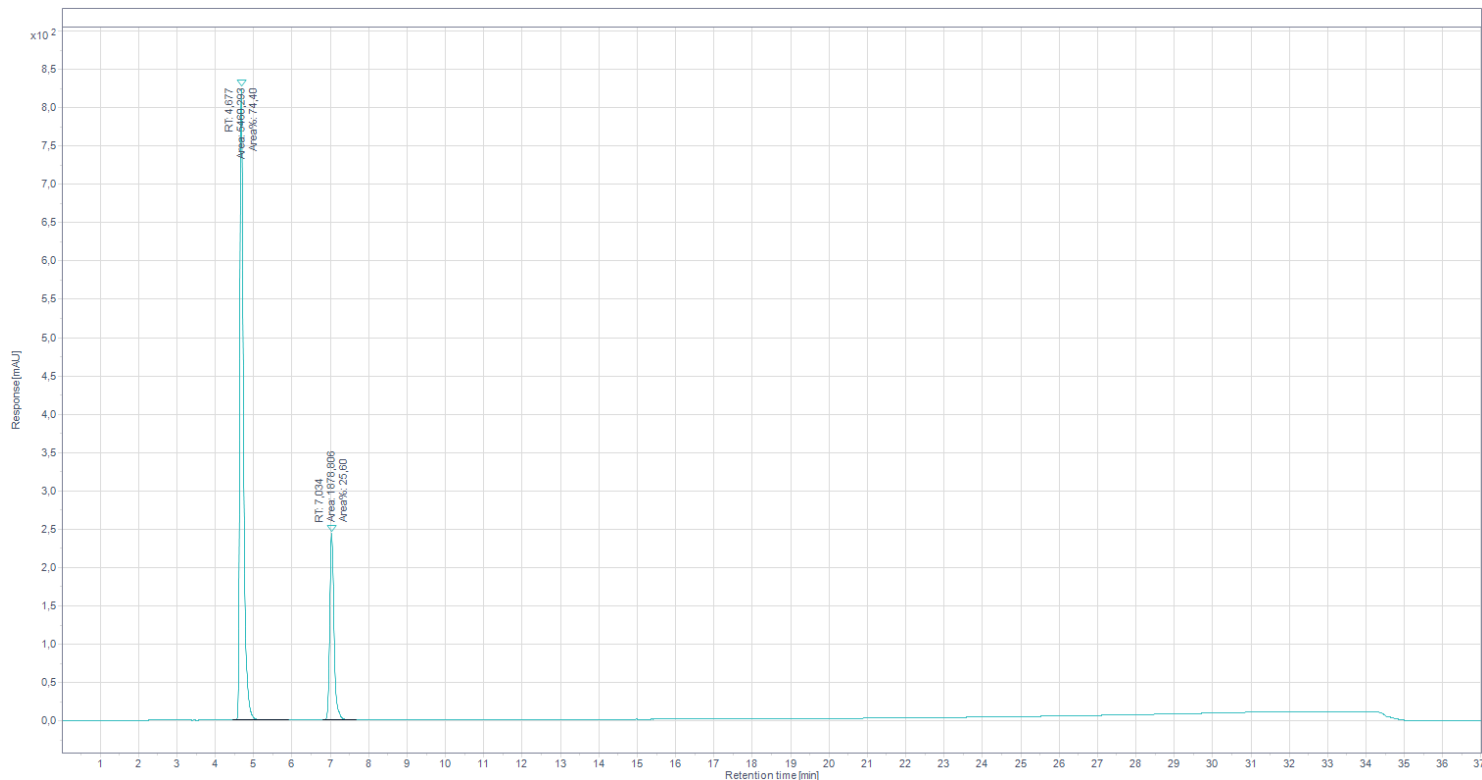

Method

Time	H2O + 0,1% TFA	ACN + 0,1% TFA
0min	60%	40%
2min	60%	40%
12min	5%	95%
15min	5%	95%

Password

D63NFU74

Analiza Białek sp. z o.o.

KRS: 0000966168

NIP: 8992922970

REGON: 521751875

Headquarters and Laboratory:

ul. Duńska 9, 54-427 Wrocław
Budynek DELTA

Lab technician:

Analysis and research specialist
mgr Arkadiusz Zajac